

basic

fixed laboratory system

The A.T. Villa BASIC system is fixed casework at its finest. Designed to outfit both small support spaces or entire laboratories, BASIC can be configured to create a custom environment with a cost conscious product.

BASIC base cabinet offerings are available in heights of 30", 34" and 36", six widths and with a myriad of door drawer configurations. A full offering of 84" tall floor cabinets and wall hung cabinets and wall shelves ensure that any space can be transformed into a functional laboratory with complementing product and aesthetic.



CORE SUPPORT



Carefully designed to integrate with our FORTE, ERGOLAB and BRAVO adaptable benching solutions, BASIC provides a cost effective, proven solution that completes the A.T. Villa product family

- Full SEFA 8 Compliance
- Removable back panels for easy service access
- Laminate, Epoxy, Phenolic and Stainless Steel work top options
- Full compatibility with A.T. Villa's CORE SUPPORT Service Distribution System




basic

fixed laboratory system



BASIC island and peninsula configurations that utilize A.T. Villa's CORE SUPPORT transform themselves into adaptable and re-configurable benching solutions. CORE SUPPORT also allows the addition of adjustable overhead shelving and storage options at any time.

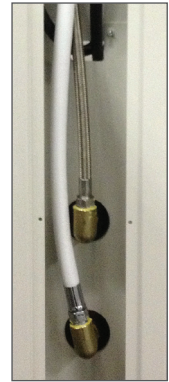

CORE SUPPORT can also be used to neatly route and organize laboratory services such as power, data, gases and plumbing. CORE SUPPORT can even be used as low profile service chase when extended to the ceiling grid.



CORE SUPPORT

The A.T. Villa BASIC and CORE SUPPORT offerings combine to provide the most cost effective and adaptable fixed benching solution available.



Widths Depths Heights

15", 18", 24", 30", 36", 42", 48", 58" end bench
14", 24", 30", 36"
30", 34" [ADA], 36" Finished Worktop Height
24", 31", 36", 48" Wall Cabinets
84" Tall Floor Cases

Finishes

Painted steel, 304SS and 316SS

Construction

Steel Box
Steel Inset
Steel and wood flush overlay

Services

UL Listed Raceways, UL Listed Power Strips, CSA Certified
Inert Gases & CSA Certified Burning Gases



www.atvilla.com
sales@atvilla.com
p 800.554.9259
f 978.582.7488

Contact us today to learn how world-class products, innovative global - design experience and a passion for excellence result in truly superior research environments.