

KEY ACCESSORIES

STABILITY / ENVIRONMENTAL CHAMBERS

Models: 7000-10/25/33/50/50/75



One-Handed Door Latch

Secure yet easy chamber access. Pull handle to open, door latches automatically. Available with or without integral lock, for 25 & 33 ft³ units.



Column Air Dryer

Permits chamber operation down to very low RH setpoints. Requires external compressed dry air source. Noise suppression kit also available for in-lab use.




Extended Temperature Range

Extends unit operation down to 0°C, includes defrost cycle for 0-4.5°C



Condensate Recirculator

Provides a source of Type II equipment water and eliminates the need for a drain line. Includes sediment & DI filters, and a 7-gallon reservoir. Compact design fits between chambers and rolls out for easy access.

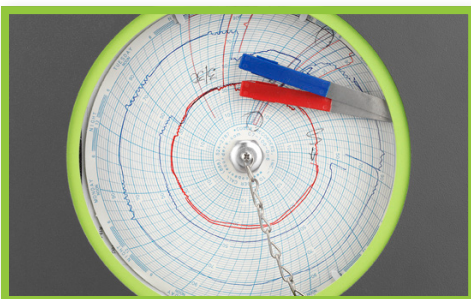

Outputs / Recorders

Circular temp/RH chart recorders in 6", 10" and 12" sizes, plus analog outputs and digital USB datalogging.