


SMARTBOX®

EA-32000 SMARTBOX® Lab Control Unit

This economical unit is designed for small scale use by individual investigators in their laboratories. The Controller automates on/off gas flow at the push of a button, delivering CO₂ to an individual cage or Euthanex Lid. The EA-32000 has a single gas output and two programmable cycles that include charge and dwell times. The run and ready lights visually indicate operation and completion.

Operates at a fixed flow rate based on the external CO₂ source.

Compliant with AVMA Guidelines

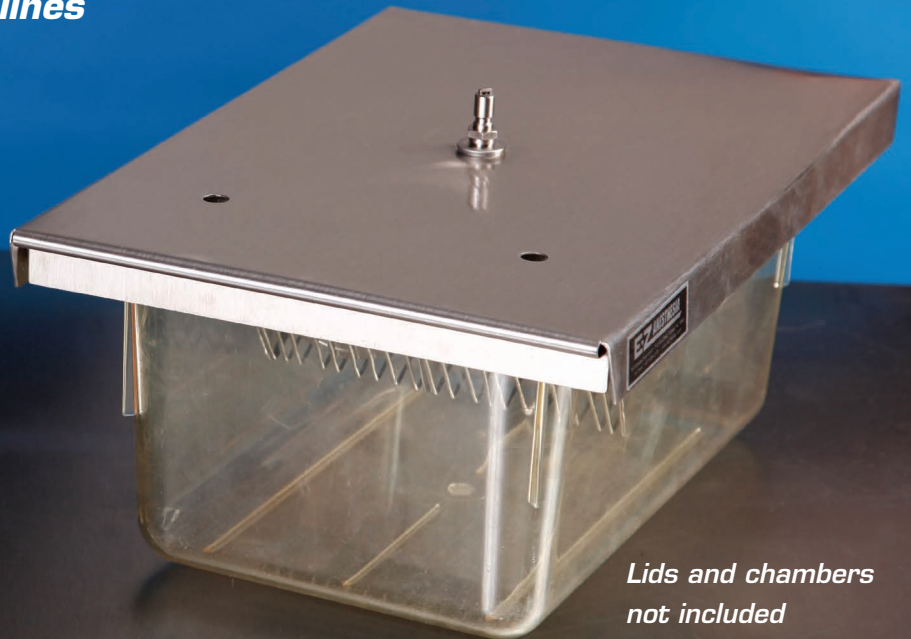


Lids and chambers not included

EA-32000 includes: SMARTBOX® Lab Control Unit and one EA-1130 Lab Unit Hose


Lab Control Unit Features:

- Treats animals at individual cage level
- Designed for investigators that need to retrieve animals quickly
- Utilizing the full range of Euthanex Lids or existing induction chamber
- Minimizes operator error and excessive gas use
- Automatic shut-off saves gas and protects personnel
- Double throw switch allows for two different treatment times
- Lights indicate cycle completion
- Easily attaches to wall

Accessories

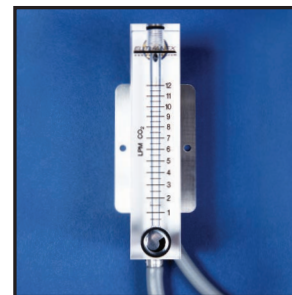

512-A1-BR

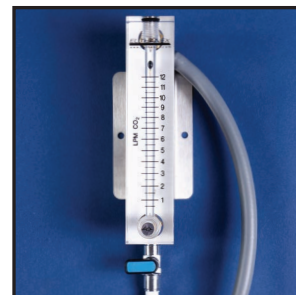

512-T1-BR




EP-PS50




EP-1305


SMARTBOX®

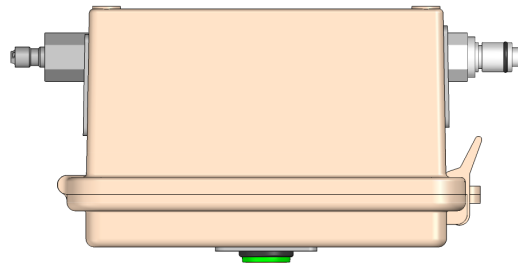


TOP VIEW

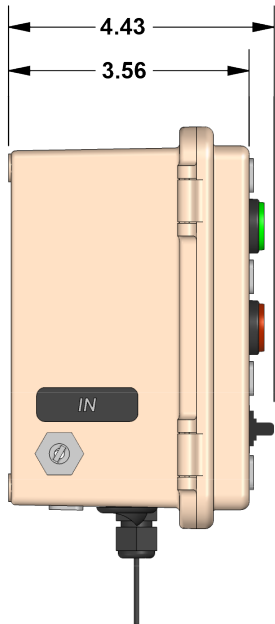


LEFT SIDE VIEW

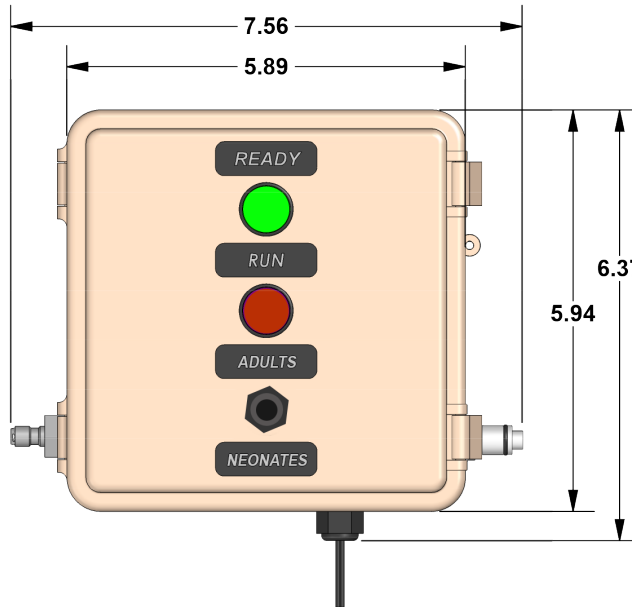

FRONT VIEW

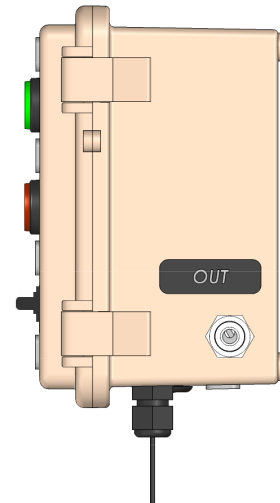

RIGHT SIDE VIEW

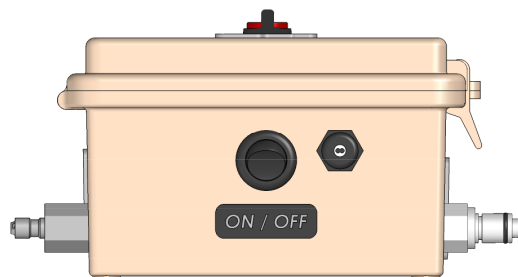

BOTTOM VIEW

SMART BOX LAB CONTROL UNIT
EA-32000

E-Z SYSTEMS
www.euthanex.com
1-877-559-0159

NOTE:
REQUIRES 110 VAC.