

EA-34004C SMARTBOX® Lab Control Unit

The Controller automates on/off gas flow at the push of a button, delivering CO₂ to individual cages or Euthanex Lids. The EA-34004C has four separate gas outputs allowing up to four cages to be euthanized individually or staggered. This feature increases the user efficiency when individually euthanizing a large number of animals. The EA-34004C has four separate switches that offer independent operation of each output. The run and ready lights easily indicate outputs in operation.

Operates at a fixed flow rate based on the external CO₂ source.

Compliant with AVMA Guidelines



Lids and chambers not included

EA-34004C includes: SMARTBOX® Lab Control Unit and four EA-1130 Lab Unit Hose

Lab Control Unit Features:

- Treats animals at individual cage level
- Designed for investigators that need to retrieve animals quickly
- Utilizing the full range of Euthanex Lids or existing induction chambers
- Increase efficiency by staggering cage start times
- Minimizes operator error and excessive gas use
- Automatic shut-off saves gas and protects personnel
- Double throw switches start and stop each output
- Lights indicate cycle completion
- Easily attaches to wall

Accessories

Lid with Flow Restrictor




EP-IL60



EP-1150-P



TOP VIEW

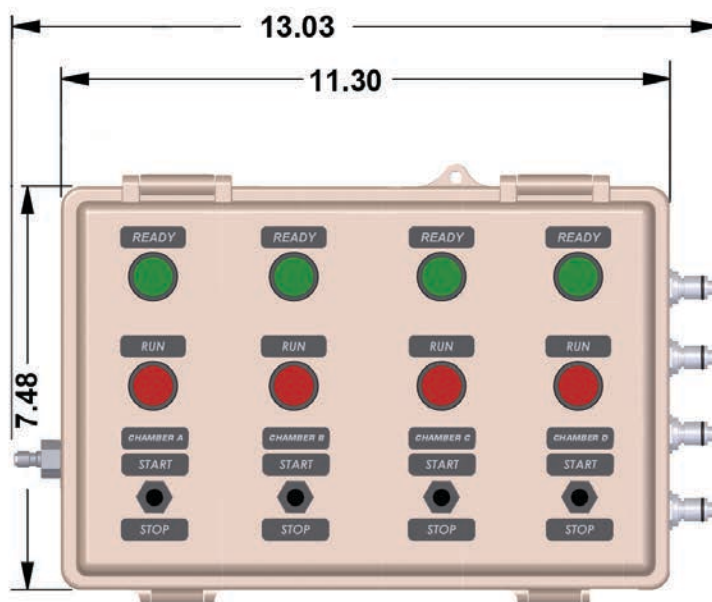

FRONT VIEW

**SMART BOX LAB CONTROL UNIT
EA-34000C**