


ALLERGARD^{ES}

NU-608 CLASS I ANIMAL BEDDING DISPOSAL CABINET

NEW! Model NU-608

Animal Bedding Disposal Cabinet

- Personnel & Environmental Protection
- Maintains an Active Air Barrier with a Constant 105 FPM [0.53 m/s] Inflow Air Velocity
- Easy-to-Clean Stainless Steel Work Surface and Side Panels
- Pneumatic Chute Lift for Easy Bag Change Access
- Stainless Steel Bang-Bar Over the Trash Chute
- Deodorizing Charcoal Exhaust Filter
- Heavy-Duty Reinforced Casters

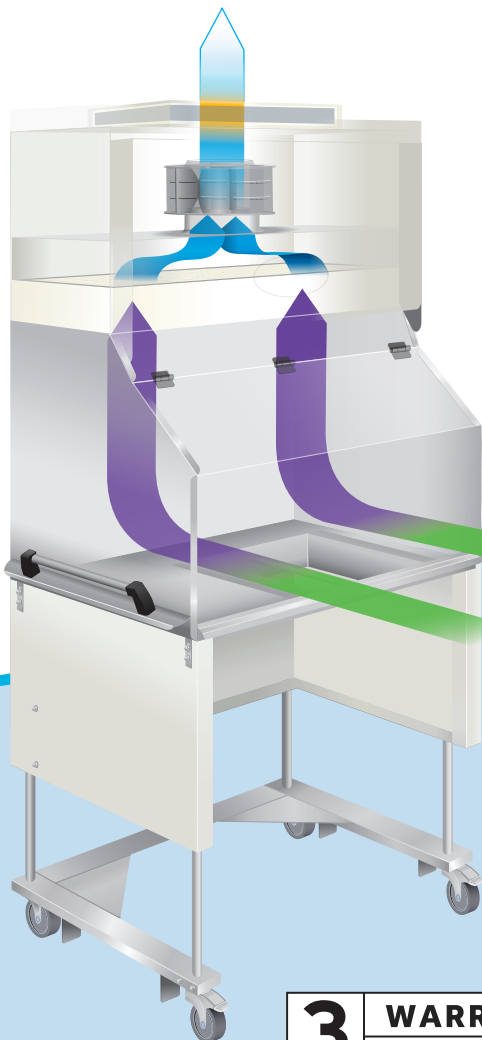

personnel &
environmental
protection

The AllerGard[®] ES (Energy Saver) NU-608 Class I Animal Bedding Disposal Cabinet minimizes the risk of exposure to particulates and allergens when performing the cage dumping process; the NU-608 makes it convenient, safe, and efficient.

3 YEAR WARRANTY
PARTS + LABOR
United States & Canada

4 YEAR WARRANTY
PARTS
Worldwide




NuAire, Inc. | 2100 Fernbrook Lane | Plymouth, MN 55447 | U.S.A. | www.nuaire.com

20-1584-F-G-EN-1-0221 © Copyright 2021. NuAire, Inc. All Rights Reserved.

ALLEREGARD^{ES}

NU-608 CLASS I ANIMAL BEDDING DISPOSAL CABINET

AeroMax™ Electronic Control System

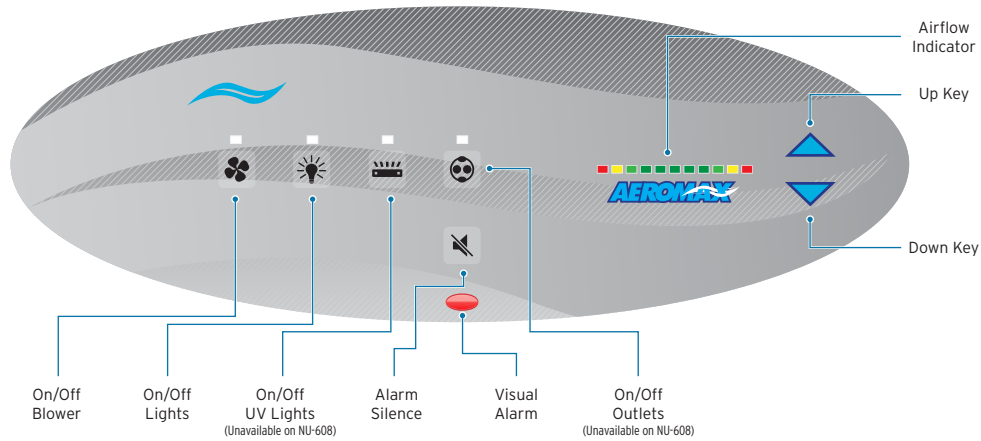

NU-608 FEATURES & SPECIFICATIONS

Class I Animal Bedding Disposal Cabinet

Model #	NU-608-400
Cabinet Style	Animal Waste Station (Console)
Protection	Personnel, Environment
Nominal Width	4 ft [1.2 m]
Workspace Dimensions (W x D x H)	46.25 x 29.75 x 28 in [1174 x 756 x 711 mm]
Exterior Dimensions (W x D x H)	52.75 x 33.4 x (84.75 - 96.75) in [1340 x 848 x (2153 - 2457) mm]
Net Weight	400 lbs [181 kg]
Certifications	UL listed, ADA Compliant
Access Opening	14 in [356 mm]
Electrical Requirements	115 VAC / 60 Hz
Air Barrier	105 fpm [0.53 m/s] Inflow
Filter Type	6 in [152 mm] Single HEPA 99.99% @ 0.3 micron Efficiency
Exhaust Module	C-25 Charcoal Exhaust Filter Module
Control System	AeroMax™ w/ PressureFlow™ Monitor
Lighting	LED (50-150 fc [538 - 1600 lux])
Noise Level	Up to 65 dba
Work Zone Construction	Single Rectangular Disposal Chute w/ Stainless Steel Bang-Bar
Base Stand Options	Telescoping Base Stand w/ Casters
	Motorized Base Stand w/ Casters
	Standard Feature
	Optional Feature



To Learn More or Request a Quote:

Visit: www2.nuaire.com/608-flyer

or Call: 1-763-553-1270

Toll Free: 1-800-328-3352

