

EA-33000TS SMARTBOX® Prodigy

The SMARTBOX Prodigy allows administrators to pre-program an unlimited number of cage sizes with designated flow rates. After initial cage setup, the user simply selects the cage icon and hits start. The Prodigy uses the cage parameters and flow rates to automatically run the cycle without further input from the user. Each user can only see the cages authorized by the administrator. Treat up to 3 cages at one time. The Touch Screen display shows the cycle stage, cycle times, and flow rate.

**Compliant with
AVMA Guidelines**


Height - 3.5"

Width - 8"


Length - 12.75"

Electrical requirement:
100-240 Volt 2amp*EA-33000TS includes:**Unit and three EA-1130 Lab Unit Hoses***PROGRAM OUT THE RISK****SMARTBOX Prodigy Features:**

- Touch screen control displays flow, cage type, stage, and times
- Automatically calculates your cage type/total volume
- Automatically controls flow rates
- Administrator can set cage flow rate and times for each cage type
- Treats animals at individual cage level
- Euthanize up to 3 cages at a time
- Secondary charges and dwells for neonates
- Unlimited amount of users with password protection
- Allows user-level permissions
- Designed for investigators that need to retrieve animals quickly
- Includes magnetic feet for easy mounting to metal safety cabinets

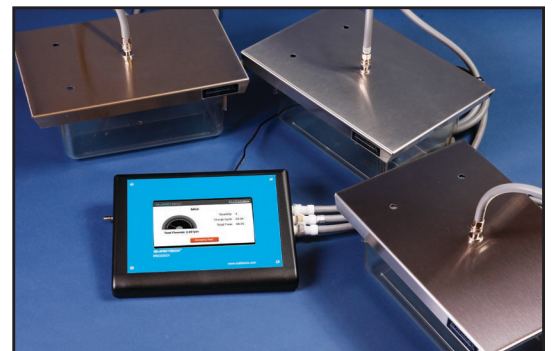

Accessories

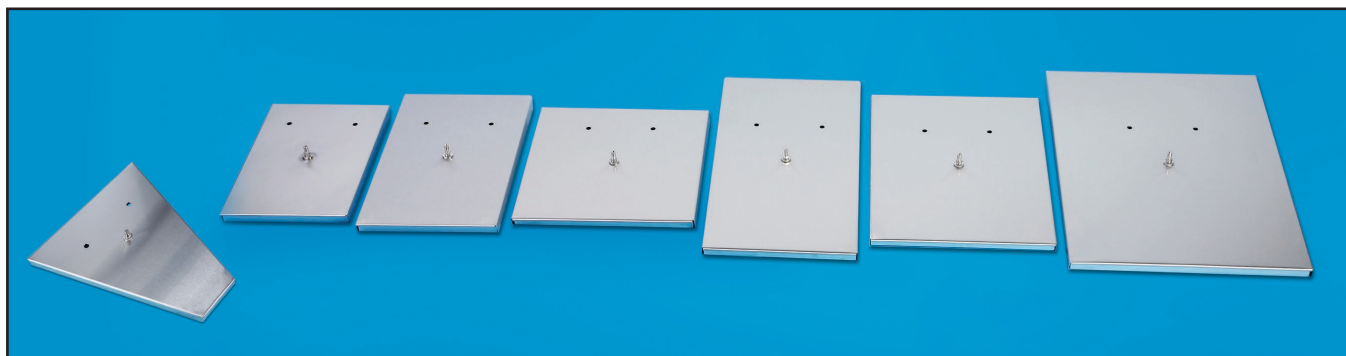

Euthanex Lids



EP-1150-P



EP-IL60




EA-550



EP-1130



Sure-Seal Induction Chambers

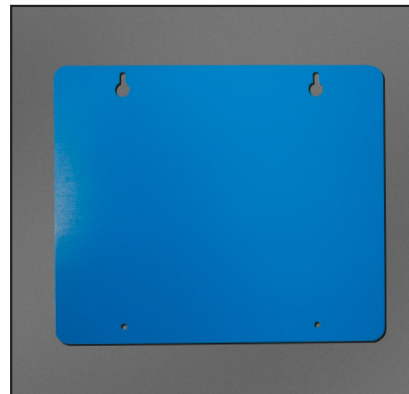

EA-1032 Metal Wall Plate
for Easy Mounting