

# PMT SCIENTIFIC

Custom  
Design

Technology

Fabrication

*Pathology Equipment  
Built Your Way*



# PMT Prodigy™

## STANDARD FEATURES

- Elevating or non-elevating styles available
- Construction is all type 304 stainless steel with #4 surface finish
- 60" long x 32" wide
- Elevating option height from 71" – 83"
- Sink – 8" deep stainless steel
- Hot/Cold water mixing faucet with swing spout and wrist blade handles.
- Hand spray – includes 5 feet of hose
- Back draft exhaust ventilation
- Full width integrated shelf located directly above exhaust grill, 8.5" deep
- Full width back splash area 52" x 22" tall for holding accessories
- (3) Large LED down lights for image quality optimization, including dimmer
- Imaging Camera Recessed Port with internal mounting hardware
- (2) Duplex GFCI protected electrical receptacles
- Ruler located on edge of counter (inch/centimeter)
- Stainless steel adjustable shelf
- Dispenser for C-fold paper towels or gloves
- Magnetic instrument tool bar
- Comes with dissection cutting board
- Unit is fully internally plumbed and wired for easy single point connection

## OPTIONAL FEATURES

- Disposal – 1 HP
- Formalin management systems
- Foot pedal activated water controls
- Removable writing ledge
- Safety splash shield
- Stainless steel pegboard
- Additional work surface lighting
- Ventilated trash receptacle
- Digital Imaging Accessory Package including computer monitor/keyboard/CPU holder and digital imaging integration
- Also available in 48" and 72" widths – including customizable for any size required!



# PMT PRODIGY FEATURES



## SEAMLESS WORK SURFACE

Keep your work surface clean and sanitary with a seamless, continuous work surface



## TOUCH SCREEN CONTROL CENTER

Monitor your grossing station's operation with work station feedback data




## ERGONOMIC & VENTILATED FORMALIN DISPENSING

Forma Fill™ Technology allows easy formalin tank fill without back strain. Formalin fumes are evacuated away with tank ventilation



## BACK DRAFT VENTILATION

Large unobstructed perforated stainless steel grill centered to allow maximum airflow efficiency and workspace

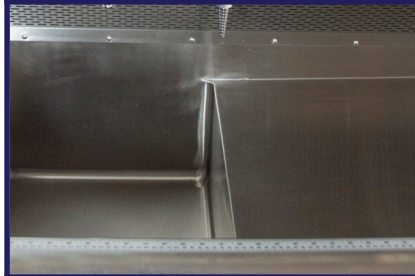

## HEAVY DUTY COMMERCIAL GRADE PLUMBING FIXTURES

All plumbing fixtures heavy brass double chrome plated for longevity



## ALL ELECTRICAL COMPONENTS CE/UL CERTIFIED

All electrical components certified to meet industrial specifications

## PMT Scientific's Safe and Innovative Formalin Management System:

PMT Scientific's latest innovation compatible with the PMT Scientific Prodigy


### Forma-Fill™ Formalin Management for collection and dispensing of formalin

- Eliminates formalin bubbling and splatter to deliver steady and even flow
- Provides ergonomic formalin transfer from carboy to ventilated tank, avoiding back strain
- High-capacity formalin tank significantly reduces time spent replenishing via carboy
- Formalin dispensing fixture is high quality and built to last
- Formalin is continuously ventilated, eliminating formalin transfer fumes
- Complete formalin level monitoring and indicators for system status at a glance



## PMT PATH STATION CUSTOMERS INCLUDE:

- Anatomic Pathology Departments
- Histology Departments
- Veterinary Pathology Facilities
- Private Pathology Labs
- Research Facilities
- Global Customer Base

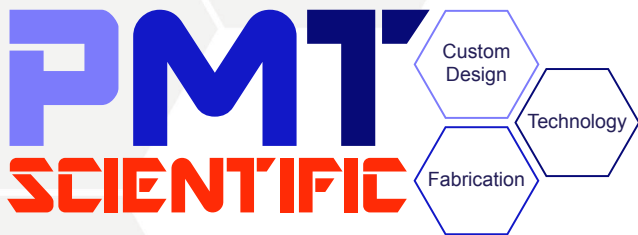

 [www.PMTScientific.com](http://www.PMTScientific.com)

 [info@pmtscientific.com](mailto:info@pmtscientific.com)

 (888) - 519 - 2286